

ahn *myDPette*

Bottle Top Dispenser

Unique piston mechanism - piston and cylinder can be cleaned without affecting the calibration


Dual graduation marks - volume can be read from both sides

Springless valve design - leads to smooth functioning and high chemical resistance

Eco-friendly - while discharging the tubes in the re-circulation mode, in fast priming, the reagent is redirected back into the bottle, eliminating any waste of precious reagents and ensuring safety and savings

Easy handling - discharge tube with 360° rotation, makes the bottle label visible all the time / telescopic filling tube compatible with different sizes of reagent bottles (bottle-height range: 170 mm to 300 mm) / ergonomic grip

Easy to calibrate - complies with ISO 8655-5 standard

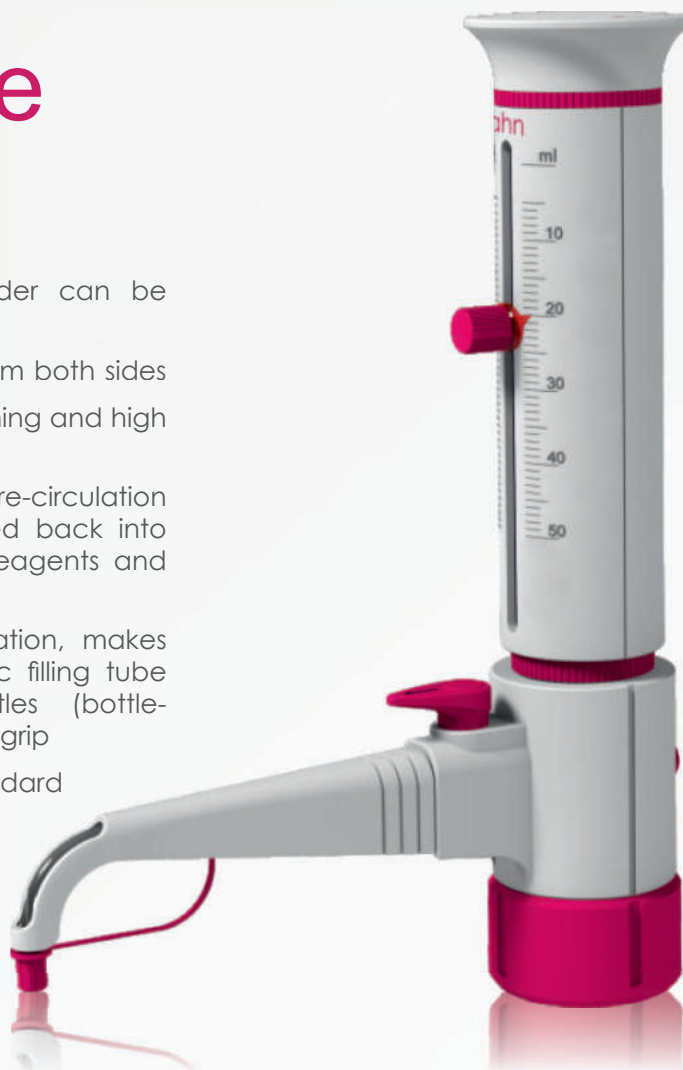

Ordering information


Description	Sales Unit	Cat. No.
myDPette Bottle Top Dispenser 0.25-2.5 mL	carton / 1 pc.	8-700-00-0
myDPette Bottle Top Dispenser 0.5-5.0 mL	carton / 1 pc.	8-701-00-0
myDPette Bottle Top Dispenser 1.0-10.0 mL	carton / 1 pc.	8-702-00-0
myDPette Bottle Top Dispenser 2.5-25.0 mL	carton / 1 pc.	8-703-00-0
myDPette Bottle Top Dispenser 5.0-50.0 mL	carton / 1 pc.	8-704-00-0
myDPette Bottle Top Dispenser 10.0-100.0 mL	carton / 1 pc.	8-705-00-0

2
YEARS
WARRANTY

autoclavable
121°C
15 min

NEW
PRODUCT

Delivery Package
1 pc. Bottle top dispenser
4 pcs. Adapters (28, 38, 40 and 45 mm screw cap)
1 pc. Telescopic filling tube
1 pc. Re-circulation tube
1 pc. Calibration tool
1 pc. User manual
1 pc. Certificate of conformity
1 pc. Calibration report
1 pc. Warranty card

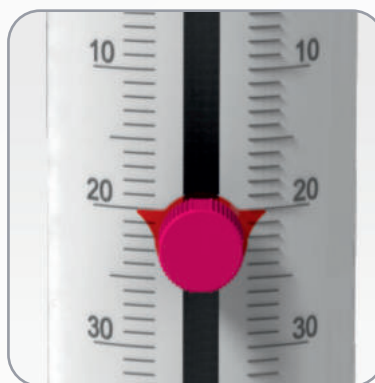

Double scale



Recirculation valve

